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Multi-Functional Dual-Output Massager


Multi-Functional Dual-Output Massager

Multi-Functional Dual-Output Massager

Multi-Functional Dual-Output Massager


Multi-Functional Dual-Output Massager

Multi-Functional Dual-Output Massager

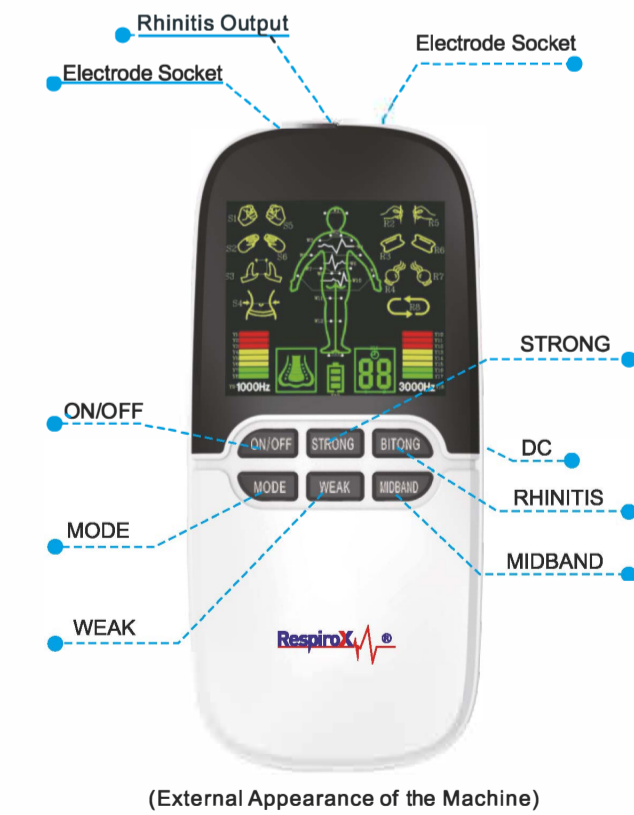(External Appearance of the Machine)

(Human Acupuncture Diagram)

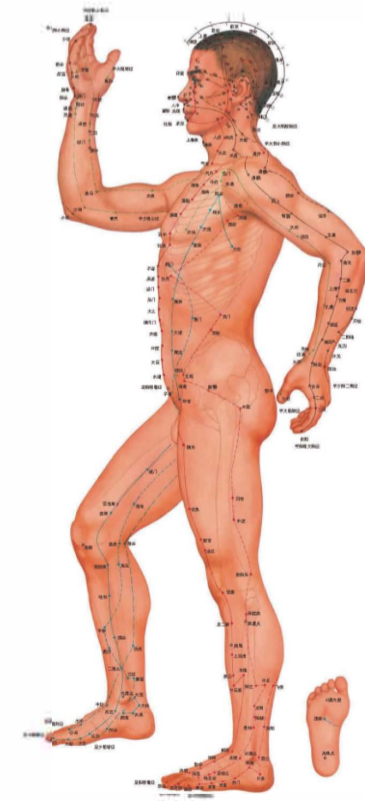Correct Way of Using:

1. Preparations
- Prepare Electrodes: Remove the patches gently, and post them on the winder when you finish using.
- Electrode Position: Wipe the treated area with a wet towel to eliminate oil or perspiration stains. Then paste up the electrodes directly on the treated areas. Bear in mind that the two electrodes should be placed on body surface at the same time otherwise the treatment cannot be carried out. Symmetrical placement of electrodes is not required.

2. Detailed Operations :
- Turn on The Power: Press the button so that the LCD is powered on and the treatment mode indicating arrow points to store mode. At the moment, the intensity is "0" and treatment time is set at 15 minutes.
- Select Treatment Mode: Press function selection button to select a treatment mode from store, massaging, manipulation, slimming, acupuncture, scraping, cupping, immune regulation and rhinitis.
- Adjust Intensity: Press the intensity button and switch the intensity to an appropriate level.

3. Precautions:

If you need to use the massager again, please press the power button and repeat the procedures above. If you change treatment mode during treatment, you have to change the electrode positions, remember to turn off the device first. If the treatment is done, please press the power button to power off the massager first before removing the electrodes. When removing electrodes, please pull the plug instead of the wire. Unavailable for pregnant women.

4. Methods for Storage and Maintenance:

- Machine cleaning: Wipe the machine with a soft cloth. Do not use diluents, gas, or volatile liquid.
- Electrode cleaning: Wipe up the stains on the electrodes with a wet cloth, or rinse them with clean water and dry them immediately before using them again.
- Storage: Keep the massager out of reach of children. Keep it away from wet places, high temperatures, direct strong sunlight, and/or water.

Configuration:

Accessories of the machine: one set of machine, two electrode wires with four electrodes, two pairs of patches, USB Line, RHINITIS Line and one instruction book.
* Electrode patches are consumables whose service life is 20-25 times. Please Prepare some extra ones.

Product Specifications:

Name	Multi-Functional Dual-Output Massager
Model	RTC-01
Power Supply	DC 4.5V
Power	17 MA (resistance load 1K)
Battery Life	8 consecutive hours for normal use
Dimensions	Length: 134mm Width:64mm, Thickness:20mm

Accessories and Maintenance

1. If you need to maintain your machine or purchase accessories, please contact the after-sale services of the manufacturer or its distributor.
2. The company warrants providing free repair services for your machine within one year ever since the purchase date in the event that any fault occurs when it is used in accordance with the instruction book (excluding artificial damages).
3. For machines exceeding the one-year warranty, the company charges only for the changed components without maintenance cost.
4. Please send the malfunctioning device to the after-sale service department or its distributor. The company will handle it immediately upon receiving it.

Shenzhen AXD Electronic Co., Ltd.
Building A2, Fuyong Zhihui Garden, Baoan District, Shenzhen City, Guangdong Province, People's Republic of China

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Date of Issue: 09.12.2022

Multi-Functional Dual-Output Massager
Model: RTC-PRO

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